

**EASY-TO-USE COMPOSITES AND BONDING AGENTS
DESIGNED FOR OPTIMAL RESTORATION RESULTS**



Natural Elegance®
Premium Flow Nano Composite



Natural Elegance®
Premium Composite



Maxima® Natural Elegance®
Packable Bulk Fill



Natural Elegance®
Flowable Bulk Fill



Natural Elegance®
Total Etch Bond



Natural Elegance®
Self Etch Plus Bond



Natural Elegance®
Universal Bond

For our **HENRY SCHEIN BRAND PRODUCTS CATALOG**, as well as additional product information, visit
WWW.HENRYSCHEINDENTAL.COM/HSBRAND.



**VALUE AND
QUALITY**



**PRODUCT
SELECTION**



**RELIABLE AND
TRUSTWORTHY**

NATURAL *Elegance*® PREMIUM FLOW NANO COMPOSITE

NATURAL ELEGANCE® PREMIUM FLOW NANO COMPOSITE is the reliable solution for smaller, hard-to-reach places and a range of common restorations. Indications for use include:

- Class III, IV, and V restorations
- Pit and fissure sealants
- Enamel defects
- Porcelain and/or veneer restorations
- Restorations for minimally invasive cavities



ADVANTAGES

- Optimized flowability make it ideal as a liner/base.
- Formulated with nanoparticles to provide strength, polish-ability, and wear resistance.
- Fluoride-releasing to prevent caries.
- Radiopacity that allows restoration visibility on x-rays.

SPECIFICATIONS

- Compressive Strength (MPa): 370
- Flexural Strength (MPa): 134
- Light Intensity for Curing (mW/cm²): 1400
- Wavelength (nm) for Curing: 450-470
- Radiopacity: 200%
- Depth of Cure (mm) 5.0
- Curing Time (sec.): 30

Premium Flow Nano Composite Syringe Kit

Contains: 4 X 2 Gm syringes (1 of each enamel shade: A1, A2, B1 & C2)
& 12 X 19 Ga syringe tips
(570-2454)

Premium Flow Nano Composite Syringe Refill

Contains: 4 X 2 Gm syringes

Enamel Shades

A1	(570-2455)	A3.5	(570-2458)
A2	(570-2456)	B1	(570-2459)
A3	(570-2457)	C2	(570-2460)

For more information, please see product Instructions for Use.

FOR OUR HENRY SCHEIN BRAND PRODUCTS CATALOG, AS WELL AS ADDITIONAL PRODUCT INFORMATION, VISIT WWW.HENRYSCHEINDENTAL.COM/HSBRAND.



**VALUE AND
QUALITY**



**EXTENSIVE PRODUCT
SELECTION**



**RELIABLE AND
TRUSTWORTHY**

NATURAL *Elegance*® PREMIUM COMPOSITE

NATURAL ELEGANCE® PREMIUM COMPOSITE delivers the ultimate in restoration quality. It contains a unique combination of photo initiators shown to provide total shade predictability and a long-lasting, smooth surface with high shine for a premier experience for both you and your patients. This universal composite is appropriate for both posterior and anterior fillings.




ADVANTAGES

- Unique mix of photo initiators
- New multi-hybrid filler and innovative photo initiator technology
- No significant color change before and after light curing
- Better mechanical properties due to improved polymerization
- Spheroidal Zircon Silicate filler gives more brightness and wear resistance
- Easy to polish with a high shine longevity result
- BPA-free and excellent non-sticky handling

SPECIFICATIONS

- Filler Rate: 73.5% – 75.5% weight (53 – 55% volume)
- Particle: Zirconium Silicate
- Particle size (Nm): 200
- Flexural Strength (MPa): 120
- Radiopacity (mm): 2.31
- Shrinkage: 0.63%
- Depth of cure range (mm): 2.30 – 3.00

AFTER 50,000 SIMULATED BRUSHING CYCLES, THE SUPERFICIAL ROUGHNESS REMAINS INTACT

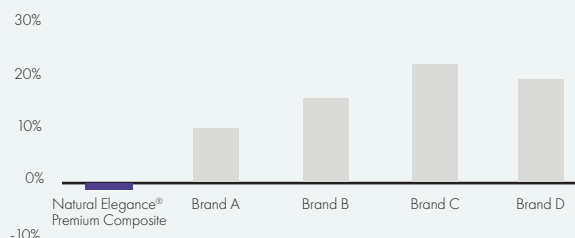

Source: Internal study

Natural Elegance® Premium Syringe

Contains: 1 x 4 Gm syringe

Enamel Shades		Enamel Shades	
A1	(570-2576)	B1	(570-2579)
A2	(570-2577)	B2	(570-2581)
A3	(570-2578)	C1	(570-2582)
A3.5	(570-2580)	C2	(570-2583)

Natural Elegance® Premium Caps

Contains: 20 x 0.2 Gm caps per box

Enamel Shades		Enamel Shades		Dentin Shades	
A1	(570-2565)	B1	(570-2569)	A1	(570-2573)
A2	(570-2566)	B2	(570-2570)	A2	(570-2574)
A3	(570-2567)	C1	(570-2571)	A3	(570-2575)
A3.5	(570-2568)	C2	(570-2572)		



For more information, please see product Instructions for Use.

FOR OUR HENRY SCHEIN BRAND PRODUCTS CATALOG, AS WELL AS ADDITIONAL PRODUCT INFORMATION, VISIT WWW.HENRYSCHEINDENTAL.COM/HSBRAND.



VALUE AND
QUALITY



PRODUCT
SELECTION



RELIABLE AND
TRUSTWORTHY

Maxima® NATURAL Elegance® BULK FILL: FLOWABLE AND PACKABLE

Henry Schein Natural Elegance® offers high-quality Flowable and Packable Bulk Fill composites so you can provide excellent aesthetic results for your patients at a great value for your practice.

- Depth of cure: 4mm • Highly radiopaque • Exposure: 10 seconds • Universal light curing • Minimal shrinkage •

Natural Elegance® Flowable Bulk Fill

- Self levelling
- Recommend for use with Natural Elegance Bond¹
- Available in syringes

Primary Indications

- Small, non-occlusal stress-bearing Class I restorations
- Base/Cavity lining under direct restorative materials in Class I and II cavities²

Technical Data

Shrinkage:	2.63%	Radiopacity:	350%
Filler weight:	75% by weight ³	Compressive strength:	357 Mpa



Maxima® Natural Elegance® Packable Bulk Fill

- Suitable for occlusal restorations and/or final layer in combination with a flowable bulk fill composite
- Recommend for use with Natural Elegance Bond¹
- Available in syringes and capsules

Primary Indications

- Class I and II posterior fillings
- Core buildup

Technical Data

Shrinkage:	1.7%	Radiopacity:	330%
Filler weight:	86% by weight ³	Compressive strength:	380 Mpa



1. Can be used with a dental bonding material of choice.

2. When used as a base/liner, Natural Elegance® Flowable Bulk Fill should be overlaid occlusally with at least a 2mm thick surface layer of a methacrylate-based universal/posterior or for instance with Maxima® Natural Elegance® Packable Bulk Fill.

3. Inorganic fillers in a methacrylate matrix (Bis-GMA, UDMA, TEGDMA).

Flowable Bulk Fill Syringe Kit

Contains: 2 x 2 Gm syringe & 20 application cannulas

Universal (900-7839)
A2 (900-7840)

Packable Bulk Fill Syringe

Contains: 1 x 5 Gm syringe

Universal (900-7691)

Packable Bulk Fill Caps

Contains: 15 x 0.25 Gm caps

Universal (900-7692)

For more information, please see product Instructions for Use.

FOR OUR HENRY SCHEIN BRAND PRODUCTS CATALOG, AS WELL AS ADDITIONAL PRODUCT INFORMATION, VISIT WWW.HENRYSCHEINDENTAL.COM/HSBRAND.



VALUE AND
QUALITY



PRODUCT
SELECTION



RELIABLE AND
TRUSTWORTHY

NATURALElegance® TOTAL ETCH BOND, SELF ETCH BOND, AND UNIVERSAL BOND

The new line of Henry Schein Natural Elegance bonding agents provides high-quality solutions at a great value.

- One component for easy application
- Less waste with drop control dispensing technology
- Compatible with all resin materials
- Universal light cured
- Solvent-free means no evaporation or drying out
- Stable at room temperature

Total Etch Bond

- **Good shear bond strength – Dentine: 22.7 MPA; Enamel: 20.9 MPA**
- Provides ultimate bonding control
- The use of etching gels is required prior to application

Self Etch Bond

- **Excellent shear bond strength – Dentine: 25.9 MPA; Enamel: 19.5 MPA**
- Saves time with fewer steps with a combined primer + bonding agent
- Minimizes post-operative sensitivity with excellent marginal sealing
- Contains 10 MDP (10 Methacryloyloxydecyl Dihydrogen Phosphate) for enhanced adhesion
- Etching of enamel and dentin is optional, not required

Universal Bond

- **Ultimate shear bond strength – Dentine: 27.5 MPA; Enamel: 22.2 MPA**
- Saves time with fewer steps with a combined primer + bonding agent
- Minimizes post-operative sensitivity with excellent marginal sealing
- Bonds direct composite resins to dentin and enamel; can also be used to prime Zirconia, Alumina, precious and non-precious metals, and porcelain
- Very low film thickness (5.5 microns)
- Etching of enamel and dentin is optional, not required



Natural Elegance® Total Etch Bond

Shear Bond Strength
Contains: 5 mL bottle
(570-2603)

Natural Elegance® Self Etch Plus Bond

Shear Bond Strength
Contains: 5 mL bottle
(570-2604)

Natural Elegance® Universal Bond Bottle

Shear Bond Strength
Contains: 5 mL bottle
(570-0603)

Natural Elegance® Universal Bond Unit Dose

Shear Bond Strength
Contains: 50 per box - 0.1 mL unit dose capsules
(570-0604)

For more information, please see product Instructions for Use.

FOR OUR HENRY SCHEIN BRAND PRODUCTS CATALOG, AS WELL AS ADDITIONAL PRODUCT INFORMATION, VISIT WWW.HENRYSCHINDENTAL.COM/HSBRAND.



VALUE AND
QUALITY



PRODUCT
SELECTION



RELIABLE AND
TRUSTWORTHY