

3M[™] Scotchbond[™] Universal Plus Adhesive

More than just a strong bond.

Universal Plus Adhesive takes universal adhesives to the next level. It keeps all the benefits of the original 3M[™] Scotchbond[™] Scotchb Universal Adhesive, the first truly universal adhesive, while adding even greater control and predictability - because a good adhesive should provide more than just a strong bond.

3M™ Scotchbond™



From excellent...

Built on the clinically proven and world's most researched universal adhesive.

> All direct and indirect indications light-curing, self-curing or dual-curing methacrylate-based composites / cements

- All etching techniques total-etch, self-etch, and selective enamel etch
- Bonds to all restorative materials without the need for separate primers
 - Available in vial and unit dose for efficient hygiene management

Virtually no post-operative sensitivity in field evaluation:

99% of dentists reported that 3M[™] Scotchbond[™] Universal Plus Adhesive shows virtually no postoperative sensitivity.*

* Field evaluation with 309 dentists and over 20,000 applications

... to exceptional

Formulation carefully optimized to achieve unique features.

First radiopague universal adhesive reducing the risk of X-ray misdiagnosis and invasive overtreatment

ersal Plus

MD

- Bonds and seals caries-affected dentin maximizing preservation of natural tooth structure, supporting minimally invasive preparation guidelines

Gold standard adhesion to all dental substrates including glass ceramics



Full dual- and self-cure compatibility

no need for an additional dual-cure activator bottle



BPA free formulation free of BPA derivatives like **BisGMA**




Fully aligned system with 3M[™] RelyX[™] Universal **Resin Cement**

excellent bond strength for virtually all resin cement indications. Adhesive cured by cement.

Added radiopacity

On average, dentists see at least one hard to interpret X-ray per day due to radiolucency under an existing restoration. This could lead to misdiagnosis and invasive overtreatment. 3M[™] Scotchbond[™] Universal Plus Adhesive is the first radiopaque universal adhesive with dentin-like radiopacity, helping to overcome this issue.

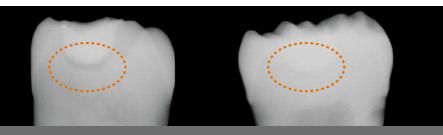

Non-Radiopaque Adhesive 3M[™] Scotchbond[™] Universal Plus Adhesive

Intentional adhesive pooling under Class I restorations – X-rays of extracted human molars, filling material 3M[™] Filtek[™] One Bulk Fill Restorative. 3M internal data, H. Loll et al., J Dent Res 99 (Spec Iss A): No. 757, 2020

Ordering information

ltem no.	Product information
3M™ Scotchbond™ Universal Plus Adhesive – Vial	
41293	Intro Kit Vial 1 Vial 5ml, 50 Disposable Brushes, 1 Etchant Syringe 3ml, 25 Dispensing Tips
41294	Refill Vial – 1 × 5ml
41295	Refill Vial – 3 × 5ml
41263	3M[™] Scotchbond[™] Universal Etchant Etching Gel 2 Syringes x 3ml, 50 Dispensing Tips

3M[™] Scotchbond[™] Universal Plus Adhesive – Unit Dose41297Intro Kit Unit Dose
50 Unit Doses, 1 Etchant Syringe 3ml, 25 Dispensing Tips41298Unit Dose 100
100 Unit Doses, 1 Etchant Syringe 3ml, 25 Dispensing Tips41299Unit Dose 200
200 Unit Doses, 1 Etchant Syringe 3ml, 25 Dispensing Tips41304Unit Dose 400
400 Unit Doses, 1 Etchant Syringe 3ml, 25 Dispensing Tips



Vial

with smooth and even surface – easy to wipe off excess

Radiolucent adhesive pooling underneath a restoration doubles the risk of restoration replacement."

Dr. John O. Burgess,

University of Alabama, USA



for efficient hygiene management – reduces the risk of cross contamination



3M.com/ScotchbondPlus

3M Oral Care 2510 Conway Avenue St. Paul, MN 55144-1000 USA Phone 1-800-634-2249 3M.com/dental

3M, Filtek, RelyX and Scotchbond are trademarks of 3M Company or 3M Deutschland GmbH. Printed in USA. © 3M 2020. All rights reserved. 70-2013-7328-2